

# 爱飞


## ASFLEX LEVEL FREE NEEDLE

聚偏二氟乙烯合成非吸收性缝合线 Polyvinylidene-Fluoride Synthetic Non-absorbable Suture



不改变缝合线的线径，  
用缝合线线径接近的缝合针与缝合线相结合的方式，  
几乎实现了消除针线结合部位的段差。

Minimized step difference of the joint between needle and suture  
by closing needle size towards suture size without changing suture size.

### 段差比较的数据 Comparing the Level at the Needle-Suture Swaged End

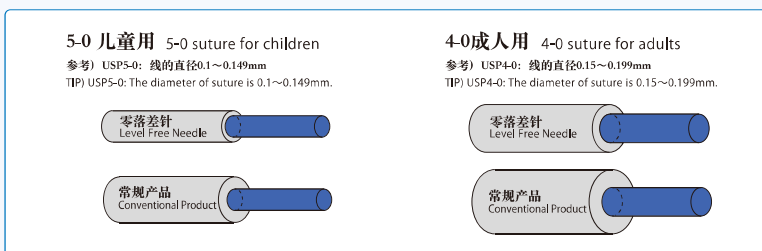


针径与线径之比 Needle-to-Suture Ratio

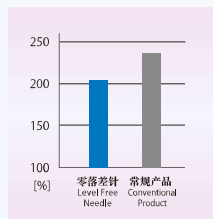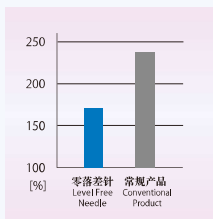

### 液体泄漏验证的数据 Data of Leak Test



5-0 儿童用

5-0 suture for children

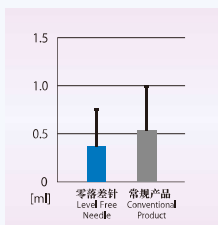

4-0 成人用

4-0 suture for adults

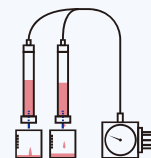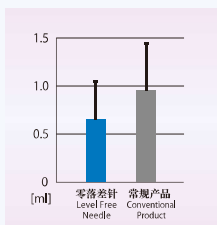








#### 泄漏实验的方法

缝合针穿过树脂板(厚度0.3mm),缝合线处于通过状态。将放入上述树脂板且装入液体的试管安装上盖子,施加一定的压力,测定一定时间内泄漏的液量。

#### Method of Leak Test

Penetrate a resin sheet (thickness 0.3mm) with the needle, with the suture placed through the sheet.  
Seal the tube with this sheet and place a liquid in it.  
Apply pressure in the tube and measure the amount of liquid that leaks out within a fixed period of time.

圆针 LEVEL FREE NEEDLE

缝合针的规格 Needle Size	数量 Quantity	颜色 Color	线长 Suture Length	线的尺寸 U.S.P.						
				7-0	6-0	5-0	4-0	3-0	2-0	1-0
<b>PLF 弱弯 3/8 Circle</b>										
 PLF-8mm 双针 Double arm	1×30袋 1pc/pkg×30	蓝 Blue	60cm		APLF716					
 PLF-9mm 双针 Double arm	4×10袋 4pc/pkg×10	蓝 Blue	60cm			APLF875S				
 PLF-11mm 双针 Double arm	1×30袋 1pc/pkg×30	蓝 Blue	75cm		APLF526					
	4×10袋 4pc/pkg×10	蓝 Blue	75cm		APLF526S	APLF705S				
 PLF-13mm 双针 Double arm	1×30袋 1pc/pkg×30	蓝 Blue	45cm			APLF505				
			60cm			APLF525				
	75cm					APLF504				
	4×10袋 4pc/pkg×10	蓝 Blue	75cm			APLF515S				
<b>RLF 强弯 1/2 Circle</b>										
 RLF-9mm 双针 Double arm	1×30袋 1pc/pkg×30	蓝 Blue	75cm		ARLF506					
 RLF-17mm 双针 Double arm	1×30袋 1pc/pkg×30	蓝 Blue	90cm					ARLF503		
 RLF-22mm 双针 Double arm	4×10袋 4pc/pkg×10	蓝 Blue	90cm				ARLF564S	ARLF533S		